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Description

BLU GEL A free is an isotonic sterile ophthalmic solution with hyaluronic acid sodium salt (FHA 1.0), and amino acids, preservative free, aimed at complementing tear secretion and restoring tear film physiology.

BLU GEL A free protects and lubricates the surface of the eye, providing relief from mechanical, environmental, and visual stress, or in the event of severe reduction in corneal-conjunctival lubrication.

Hyaluronic acid is a mucopolysaccharide physiologically present in the human eye and has the property of adhering to the surface of the eye protecting and lubricating it. Thanks to its mucoadhesive properties, it remains on the surface of the eye for a long time and stabilises the tear film.

The amino acid formulation helps lubricate the surface of the eye and improve the protection of epithelia subject to various kinds of stresses.

The BLU GEL A free bottle is fitted with an OSD (Ophthalmic Squeeze Dispenser) dropper which, through its multi-layer air filter and tip-seal valve system, ensures that the sterility of the product is maintained without preservatives.

Indications

BLU GEL A free is indicated in severe tear dysfunction syndrome for hydration and lubrication of the surface of the eye.

Dosing regimen and method of administration

Apply 1-2 drops to the conjunctival sac 3 or more times a day.


The product may be instilled into the eyes even when wearing contact lenses.

Contraindications

Known hypersensitivity to product components or other substances that are chemically closely related.

Instructions for use

- When first used, remove the plastic ring that seals the cap by pulling the side tab (1).
- Remove the cap (2).
- Turn the device over the eye and apply increasing pressure on the bottle with your index finger and thumb until the drop (3) comes out. As the solution in the bottle becomes depleted, more force may be required (4).
- Close the bottle again with the cap (5).

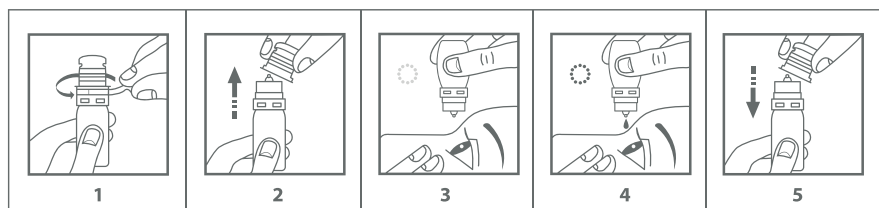Precautions and warnings

- This product is for external ophthalmic use only.
- To eliminate the risk of contamination of the solution, make sure the dropper does not come into contact with the eye or other objects;
- To avoid damage to the surface of the eye, do not allow the dropper to come into contact with the eye during instillation.
- As the solution in the bottle runs out, more force will need to be applied to create a drop.
- Close the bottle tightly after use.
- If used alongside topical ophthalmic medicines, consult your doctor.
- Do not use the product if the OSD bottle is damaged.
- Keep the product out of the sight and reach of children.

Side effects

No undesirable effects have been reported after administration of the product. In case of adverse reactions, stop using the product and consult a physician. After application of the product, there may be a slight burning of the eye and temporary visual discomfort.

Ingredients

Main component: Sodium Hyaluronate 0.30% (FHA 1.0)

Other components: Boric acid, Disodium tetraborate decahydrate, Sodium chloride, Glycine, L-Leucine, L-Lysine hydrochloride, L-Proline, Water for injection.

Storage conditions

Do not store above 25°C. Do not use the product after the expiry date indicated on the bottle label.

After first opening, the product may be used until the expiry date indicated on the packaging.

Do not discard in the environment.

Leaflet last revision date: June 2020

Caution

Consult instructions for use

Use by

Batch code

Do not use if package is damaged

Sterile processed using aseptic techniques

Upper temperature limit

Manufacturer

Do not discard in the environment after use.

Fidia Farmaceutici S.p.A.

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35031 Abano Terme (PD) Italia

Tel. +39 049 8232111

www.fidiapharma.com